ST RAMPS

Art. ST-091, ST-122, ST-152

INSTRUCTIONS FOR USE

Mobilex A/S

Gronlandsvej 5, 8660 Skanderborg,
Denmark,

www.mobilex.dk

Mobilex Sp. z o. o.

ul. Duńska 3/5, 91-204 Łódź
tel.: (0-42) 214-29-11

www.mobilex.pl

Table of Contents

Language	Page	Language	Page
 English	4	 Italiano	28
 Български	6	 Lietuvių	30
 Čeština	8	 Latviešu	32
 Dansk	10	 Norsk	34
 Deutsch	12	 Nederlands	36
 Ελληνικά	14	 Polski	38
 Español	16	 Português	40
 Eesti	18	 Română	42
 Suomi	20	 Slovenčina	44
 Français	22	 Slovenščina	46
 Hrvatski	24	 Svenska	48
 Magyar	26		

ST RAMPS Illustrations

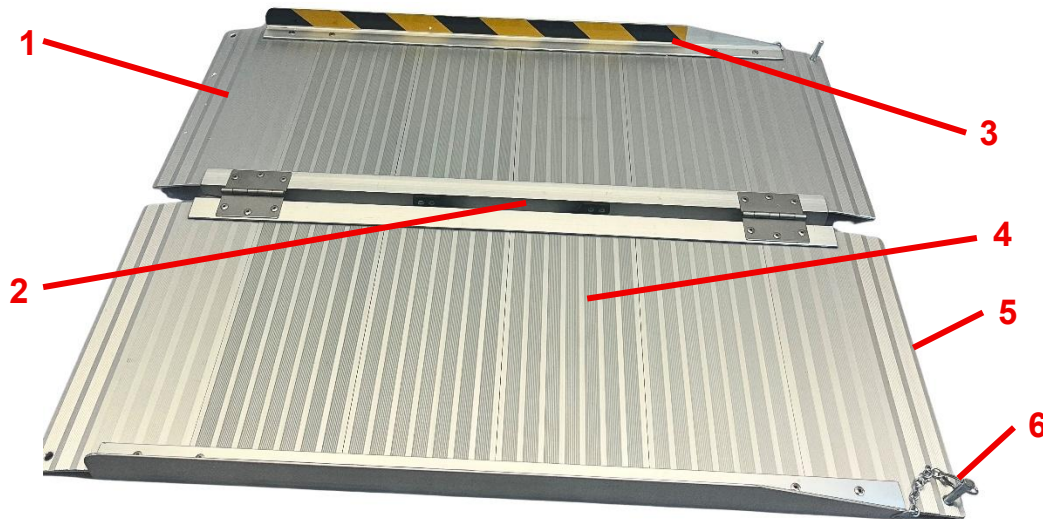
Fig. 1 – Device elements

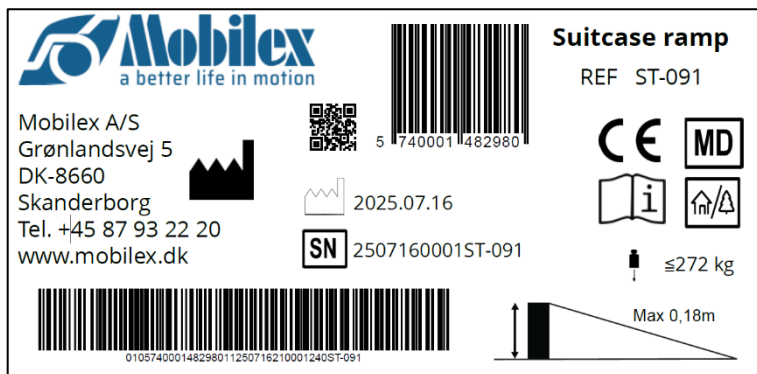
Fig. 2 – Example of the product label

	1. QR-code link to homepage		7. CE marked product
	2. Manufacturer		8. Manufacturing date
	3. Observe the user manual		9. Max load
	4. Serial no.		10. Medical Device
	5. Item Number		11. To be used indoor and outdoor
			6. Max Elevation / height

Fig. 3 – Meaning of the symbols on the product label

INSTRUCTIONS FOR USE ST Ramps - Art. ST-091 – ST-152

1. Introduction - Thank you for buying a high-quality Mobilex product. These instructions for use contain a description and important information about the use of the medical device. They are intended to help you to use the product correctly and safely. It is very important to read the instructions carefully before using the device. Pay particular attention to the safety guidelines and always follow them. Mobilex is constantly trying to improve and upgrade its products. Therefore we reserve the right to change specifications and technical data of the product without further notice. All figures referenced in this manual are on the Illustrations page (page 3).

Intended use of the product - To be used to assist people in wheelchairs or scooters by enabling the passage over elevated points.

Indications - Ramps are indicated for people with mobility limitations who due to problems with walking, standing or sitting are forced to use mobility devices such as wheelchairs, scooters, rollators etc. Ramps are indicated for elimination of architectural obstacles by smooth connection of the level differences that a wheelchair user has to overcome and thus they allow the wheelchair user to reach places normally inaccessible to him/her. The ramp is intended for use independently or with the help of an attendant. Using the ramp for purposes other than those indicated above, is prohibited.

Contraindications - Contraindications cannot be tied to a specific clinical picture. The choice of a suitable ramp type and length must always be made with the help of a specialist doctor or physiotherapist. Observe the safety information in chapter 2.

Areas of use - The product can be used in home care, health and long-term care facilities, and hospitals.

Medical incidents and unwanted events - In case of any unwanted events related to the device, please contact the dealer or manufacturer.

2. Safety requirements

- Make sure that this manual has been read by all persons using the device.
- The manufacturer assumes no responsibility for damage and/or injury resulting from failure to follow the operating instructions.
- Only use the product in a flawless condition.
- If defects or errors are detected you must immediately contact the dealer.
- Follow all manufacturer's recommendations and warnings.
- Be aware of the information on the product label.
- Only use the product for the described purpose.
- Avoid making structural changes to the device, unless you have the manufacturer's written acceptance for such modifications.
- During different adjustments which are possible on the device the user must take care not to squeeze the fingers or other body parts between two moving elements of the device.
- The device must only be used on a stable surface.
- Ramps must not be loaded with more than the "max load" described in Technical Data and on the label.
- The ramps can be used to overcome obstacles at maximum heights of: ST-091 18 cm, ST-122 24 cm, ST-152 30 cm.
- Be careful during folding or unfolding of the ramp. There is a risk of a finger getting stuck or cut between parts of the ramp.
- Check the stability and correct positioning of the ramp before each use.
- Do not use the ramp in adverse weather conditions (ice, heavy rain) that could make the surface slippery.

5. Technical Data

ST Ramps	ST-091	ST-122	ST-152
Length	91.4 cm	121.9 cm	152.4 cm
Width		73.8 cm	
Surface width		73 cm	
Height of side protect.		3.5 cm	
Max. Elevation / height	12° / 18 cm	12° / 24 cm	12° / 30 cm
Weight	8 kg	10.4 kg	13.1 kg
Max. load		272 kg	
Folded size cm (L x W x H)	101 x 43 x 10.5	131 x 43 x 11	162.5 x 43 x 11

Product label - An example of the product label is shown in Fig. 2.

Meaning of the symbols on the label - see Fig. 3, descriptions of symbols 1–11 are provided below:

- | | | | |
|-----------------------------|---------------------------|-----------------------|-----------------------------------|
| 1. QR-code link to homepage | 4. Serial no. | 7. CE marked product | 10. Medical Device |
| 2. Manufacturer | 5. Item Number | 8. Manufacturing date | 11. To be used indoor and outdoor |
| 3. Observe the user manual | 6. Max Elevation / height | 9. Max load | |

Parts supplied - Each product package should contain: 1 ramp, 1 user manual.

Check the contents of the box against the above list. If at this stage you discover that a part is missing or damaged, contact your supplier immediately.

3. Description

Product description - ST-Ramps are foldable platforms made of aluminium. They have an anti-slipping rib shaped surface.

Elements - The elements of the device are shown in Fig. 1. The descriptions of elements 1-6 are provided below:

1. Top support
2. Carrying handle (under the ramp)
3. Side protection
4. Anti-slip rib shaped surface
5. Bottom support
6. Locking mechanism (prevents the ramp from sliding away)

Function description

Unfolding and positioning the ramp

Place down the ramp and unfold it so that the carrying handle is on the underside of it.

To fold it again, simply pull one side onto the other so that the carrying handle is reachable.

Positioning the ramp

Place the top support on the elevated surface, and the bottom support on the ground level. Make sure the ramp is steady and placed on a stable surface before use.

Assembly

This device requires no assembly instructions.

Accessories - There are no accessories provided for this product.

4. Care & Maintenance - Normal dirt can be removed with standard cleaning products and a sponge or soft rag. Check the specific product information and use only commercial cleaning products which are suitable for cleaning and disinfection (do not use alcohol-containing products and no solvents or abrasives).

The device basically does not need maintenance. Despite the solid construction and the use of resistant materials the product is subject to wear. It is therefore recommended to check the product at regular intervals by a professional service. We recommend the following to be checked:

- Cleaning and disinfection of the product
- Frame/product should be checked if it has been bent, damaged, worn or corroded
- Wear parts are checked and replaced if necessary
- Checking the folding elements and the carrying handle for proper function
- Inspection of the anti-slip surface for damage or wear

Expected service life - Normally 10 years depending on use.

Storage - We recommend storage in a dry place at a temperature above 0 °C, away from intense sunlight.

Transportation - The ramp can be locked in a folded position. It should be transported in a stabilized and secured way with no risk of sudden or dangerous movements during transport.

Reuse - The product is suitable for reuse. Before a new user takes over the product, it must be checked according to point 4.

6. Most common problems and solutions

Problem	Possible cause	Solution
Ramp is not stable	The surface under the ramp is not flat	Check the position of the ramp. Put it on a flat, smooth surface
Ramp is not stable	Ramp is deployed upside down	Check if the ramp is positioned correctly
Ramp doesn't open	The latch is not released	Release the latch
Ramp does not fold	Obstacles between the folding parts	Remove dirt or foreign objects from the folding areas
Surface is slippery	Moisture, ice or contamination	Clean and dry the surface before use

7. Disposal of the product - The product cannot be disposed of with household waste but must be brought to the local recycling center. The ramp consists mainly of aluminium and can be recycled. Observe the applicable local regulations for the disposal of aluminium products.

8. Guarantee terms - Mobilex A/S offers 2 years of warranty for damages caused by errors in production or material errors. The warranty is not valid in case of not-intended repair or usage of the product. Parts subject to normal wear are not covered by the warranty, unless wear is caused by a manufacturing fault. If product changes are made without our written acceptance the CE-marking and product warranty will not be valid. In the event of product damages covered by this warranty, please notify your dealer or Mobilex A/S directly. The warranty does not cover transport costs and does not include compensation for personal injuries or non-use of the product during repair. The warranty does not cover damages due to the user manual not being followed.

9. Declaration of conformity - The device has passed the conformity assessment in accordance with EU Regulation: MDR 2017/745 regarding medical devices and is CE marked. A declaration of conformity has been issued for the device and is available at: mobilex.dk/medias/2025-12/ce_st_ramper_19935.pdf

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